

EXPERIENCE THE HEALING NATURE OF LIGHT



NEW HIGH POWERED MEDX LASER & LIGHT THERAPY TREATMENT SYSTEM

MedX® Laser and Light therapy is a proven and highly effective treatment which increases local microcirculation: accelerates tissue repair, reduces pain, decreases inflammation/ edema, and enhances immune response and angiogenesis.

The MedX Laser and Light therapy treatment system is powered by a compact MedX1100N console. The console is easy to operate with the option of using the preset guides which automatically calculate and display treatment times or the manual mode referring to the MedX Treatment Manual. The detailed MedX Treatment Manual provides accurate protocols and settings for over 65 common ailments. The console is quick and easy to set and makes preparing for treatments simple.

Three accessory ports provide the flexibility to use any combination of accessories with either continuous or pulsed waveforms. The console powers the MedX hand held low level laser and two hands-free strap-on superluminous diode (SLD) clusters. Once activated, the console displays the time and energy delivered for each port separately and automatically turns itself off once treatment is completed. Multiple sites can be treated at the same time.

MedX products are FDA cleared and Health Canada Licensed under following numbers: 15556, 76086, 69680 & 60952. Our products are produced in an ISO13485: 2003 certified manufacturing and testing facility in Canada. For a complete profile of MedX Health Corp. and our products visit www.medXhealth.com

MedX Laser and Light therapy stimulates the body's natural repair response at the cellular level to reduce pain and accelerate healing.

The MedX Laser and Light Therapy system delivers up to **2450mW** of total power



Bringing drug-free, non-invasive solutions to light

Kosim Medical Supplies


55 Valleywood Dr,
Markham, ON L3R 5L9
Tel: 1-877-50-KOSIM

Visit: www.kosim.ca or sales@kosim.ca

Kosim
ideal for health



EXPERIENCE THE HEALING NATURE OF LIGHT

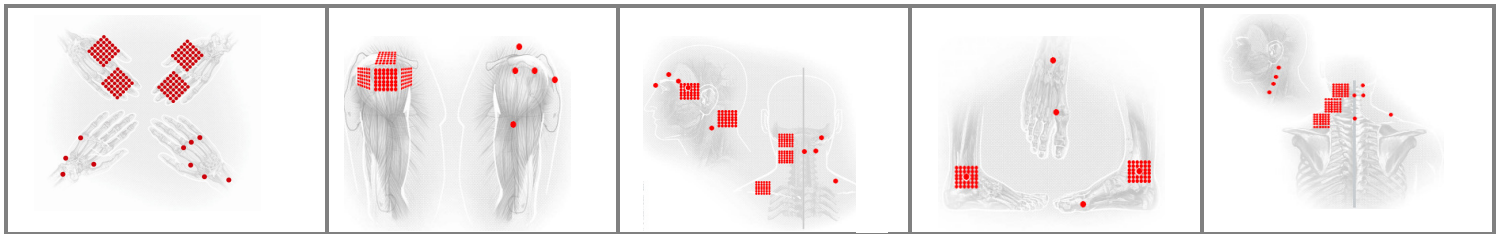

MEDX LASER & LIGHT THERAPY TREATMENT SYSTEM **INCLUDE:**

MedX1100N console:

A compact unit that fits on your treatment cart and is easy to relocate. Three accessory ports allow flexible options enabling the treatment of large areas or several sites simultaneously.

Hand-held low-level laser:

The MedX laser provides a total of 450mW of power with three 150mW GaAlAs (808nm) laser diodes and a blue visible guide LED. (Class IIIb)

SLD clusters:

The large area superluminescent diode cluster provides a total of 1000mW of energy with a combination of predominantly infrared (870nm) and visible red (633nm) diodes. Straps-on for convenient hands-free treatment.

Treatment manual:

An outline of highly effective and efficient protocols for the delivery of SLD and laser treatments for over 65 conditions.

Operating manual

Two sets of goggles

Power supply

MedX Kits: Product Description

0915	MedX 450mW Portable Laser
3116	MedX Console + 450mW Laser + 2 SLDs
2031	MedX ^{HOME} System



©MedX Health Corp. 2005



1-888-363-3112 (Cdn/US only) • 905-670-4428 • infor@medXhealth.com • www.medXhealth.com